



Testing Equipment Manufacturer
For Paint, Coating & Ink

www.bevsinfo.com

BEVS Industrial Co., Ltd. is the ISO approved laboratory testing instruments manufacturer in the field of coating, paint, resin & ink. Our quality control and testing instruments are made within the highest standards and combined with today's most advanced technology. BEVS will be the best partner with you and provide the whole solution of testing, inspection and laboratory instruments.





Application of Film

Sagging Applicator

Introduction:

The instrument is used for sag resistance evaluation to ensure a uniform film thickness by selecting suitable thickness, which is conducive to verticals, corners, holes, circular pieces of painting controlled.

Viscosity, fineness, thickness are major factors in the impact of sagging.

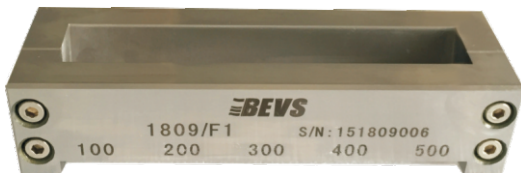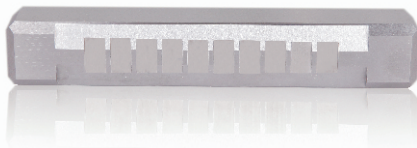

Technical Specification:

- ◆ High-grade stainless steel
- ◆ 10 different films incremented by 25µm
- ◆ Wet film width 6mm, Spacing 1.5mm

Standard:

ASTM D4400, D373

Order Information:

Order No.	Range (µm)
BEVS 1809/1 Sagging Applicator	50-75-100-125-150-175-200-225-250-275
BEVS 1809/2 Sagging Applicator	250-275-300-325-350-375-400-425-450-475
BEVS 1809/3 Sagging Applicator	450-475-500-525-550-575-600-625-650-675
BEVS 1809/4 Sagging Applicator	650-675-700-725-750-775-800-825-850-875
BEVS 1809/5 Sagging Applicator	850-875-900-925-950-975-1000-1025-1050-1075
BEVS 1809/F1 Quadruplex Sagging Applicator	100-200-300-400-500
BEVS 1809/F2 Quadruplex Sagging Applicator	600-700-800-900-1000

Other sizes are available upon request.